



# Terumo® Dental Needle

Single use



**Contamination-proof  
and clean-cutting  
for anaesthesia  
comfortable to the patient**

# Terumo® Dental Needle

## Wide selection for optimal length

TERUMO® Dental Needle is available in both short and long sizes for conduction anaesthesia and infiltration anaesthesia.

## Minimal patient discomfort

The ultra-sharp LANCET POINT minimises pain and discomfort, and eliminates the need for a surface anaesthesia prior to injection. The high-tensile TERUMO® needle is difficult to break even if the patient moves violently during injection.



## Individual containers

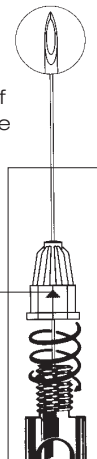
Each needle is individually packed in its own plastic sheath. The colour-coded sealing label provides instant identification of needle gauge and facilitates inventory control.

## Reliable fixing

The needle hub screws directly and securely onto the syringe. It will neither come off accidentally nor leak. Since the needle can be screwed onto the syringe without being taken out of its sheath, sterility is maintained up to the moment of injection.

### Bevel Indicator

The arrow on the hub allows you to align the bevel of the needle for easier insertion.

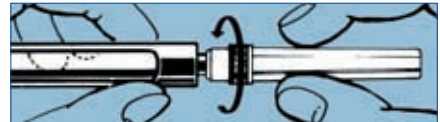


## Instructions for use

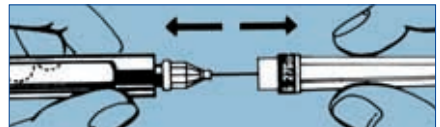
1. Twist cap to break sealing label. Remove cap.



2. Attach needle to syringe by using cartridge sheath as a wrench to screw the needle hub onto the syringe.


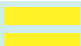



3. Immediately prior to injection, remove sheath to expose needle.



## Other features include:

- screw hub locks needle firmly onto syringe
- colour coding identifies needle gauge
- paper seal ensures sterility

Gauge	Colour code	IMPERIAL THREAD - Effective needle length							METRIC THREAD - Effective needle length				
		16 mm	19 mm	21 mm	22 mm	30 mm	35 mm	41 mm	8 mm	13 mm	16 mm	21 mm	35 mm
25G	 Red			DN*2521		DN*2530		DN*2541	DN*2508F			DN*2521B	DN*2535B
												DN*2535F	
27G	 Yellow			DN*2721	DN*2722	DN*2730	DN*2735	DN*2741			DN*2716F	DN*2721B	DN*2735B
												DN*2721F	
30G	 Green	DN*3016	DN*3019	DN*3021	DN*3022					DN*3013B	DN*3016B	DN*3021B	
										DN*3013F	DN*3016F	DN*3021F	

Needle length cartridge side is always 10.7 mm

Codes with B = needle length cartridge side is 13.2 mm  
Codes with F = needle length cartridge side is 10.7 mm

TERUMO Dental Needle is available with both Metric (blue hub) and Imperial (white hub) threads. For detailed ordering information, please contact your usual supplier or the nearest TERUMO sales office. Package: 100 pcs



Every day, as healthcare professionals, you are working to ensure life and to contribute to the well-being of your patients. Every day, Terumo Europe provides you with the clinical systems, means and solutions to keep life flowing.



We keep life flowing

# TERUMO®

**TERUMO CORPORATION**  
44-1, 2-chome  
Hatagaya  
Shibuya-ku  
Tokyo 151  
JAPAN  
Tel. +81/3 3374.8111  
Fax +81/3 3374.8196

**TERUMO EUROPE N.V.**  
Researchpark Zone 2 Haasrode  
Interleuvenlaan 40  
B-3001 Leuven  
BELGIUM  
Tel. +32/16.38.12.11  
Fax +32/16.40.02.49

**TERUMO EUROPE N.V. BENELUX SALES DIVISION**  
Researchpark Zone 2 Haasrode  
Interleuvenlaan 40  
B-3001 Leuven  
BELGIUM  
Tel. +32/16.38.14.02  
Fax +32/16.38.15.55  
THE NETHERLANDS  
Tel. +0800/023.19.38  
Fax +0800/022.34.89

**LABORATOIRES TERUMO FRANCE S.A.**  
1, Parc Ariane  
Bâtiment Uranus  
Rue Hélène Boucher  
73284 Guyancourt Cedex  
FRANCE  
Tel. +33/1 30.96.13.00  
Fax +33/1 30.43.60.85

**TERUMO (DEUTSCHLAND) GmbH**  
Lyonerstrasse 11 a  
60526 Frankfurt am Main  
GERMANY  
Tel. +49/69.66.44.20  
Fax +49/69.66.66.826

**TERUMO EUROPE SUCURSAL EN ESPAÑA**  
Avenida de Burgos 16-D  
Edificio EUROMOR  
Madrid 28036  
ESPAÑA  
Tel. +34/91 383.90.69  
+34/91 383.91.61  
Fax +34/1 767.04.60

**TERUMO EUROPE ROME BRANCH**  
Via Simone Martini 143/145  
I-00142 Roma  
ITALIA  
Tel. +39/06.519.61.420  
Fax +39/06.503.04.07

**TERUMO EUROPE N.V. Belgium filial Sverige**  
Nya Varvet  
Byggnad 90  
SE-426 71 Västra Frölunda  
SWEDEN  
Tel. +46/31.748.58.80  
Fax +46/31.748.58.90

**TERUMO U.K.**  
(Branch of Terumo Europe N.V.)  
3 Unity Grove  
Off School Lane  
Knowsley Industrial Park South  
Knowsley  
Merseyside  
L34 9GT  
UNITED KINGDOM  
Tel. +44/151.482.4800  
Fax +44/151.482.4811

®: Registered Trademark  
Published by Terumo Europe N.V.

[www.terumo.com](http://www.terumo.com)  
[www.terumo-europe.com](http://www.terumo-europe.com)